


Reflectors for R- & S-Series Probes




Probe reflectors are stainless steel in 3 configurations: L, J, U

Non-magnetic Brass & Macor Reflectors are available for MRI applications.

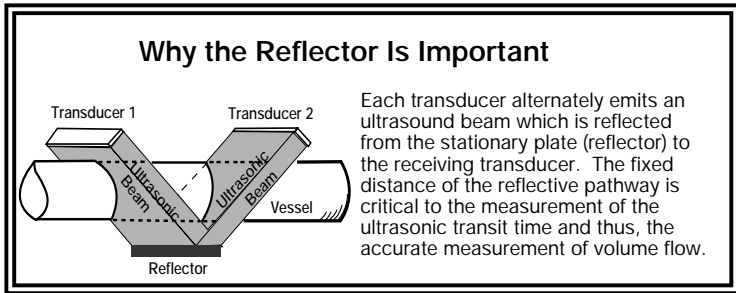

L-Reflectors


The L-shape reflector has an adjustable slide cover to fully enclose the vessel. Suture holes are provided to stabilize the probe for prolonged applications.

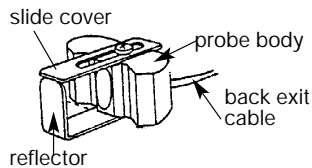



J-Reflectors

The J-shape reflector is available with or without a slide cover; useful for undersized or inaccessible vessels.

with Slide Cover

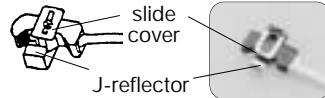


Open J reflector

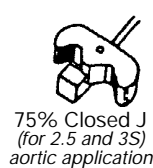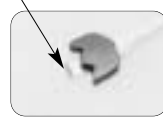

U-Reflectors

The U-shape reflector is equipped with screw mounts on each side of the probe body. To install around vessel, the reflector is removed from probe body, placed around vessel and is again secured with screws to the probe body. The reflector can be encapsulated with a silicone shield at extra cost. *(Description on opposite page.)*

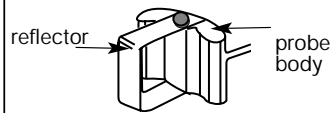

Reflector Availability

Size	L	J	U	SF	US
1R, 1.5R, 2S, 2.5S	–	a	–	–	–
2R	a	a	a	a	–
3S, 4S	a	a	a	a	–
3R, 4R	a	a	a	a	–
6S, 8S	a	a	a	a	a
6R, 8R	a	a	a	a	a
10S - 24S	a	a	a	–	a

SF - silicone flange; US - silicone encapsulated U; (a) available; (–) not available